

Anti-Human ZCCHC7, monoclonal

Alternate Names: Zinc finger CCHC domain-containing protein 7

Cat. No. 14-020

FOR RESEARCH USE ONLY

NOT FOR USE IN HUMANS



PRODUCT DESCRIPTION

UniProt Summary

UniProt

Primary Accession: [Q8N3Z6](#)

Component of a nucleolar TRAMP-like complex, an ATP-dependent exosome regulatory complex consisting of an helicase (SKIV2L2/MTR4), an oligadenylate polymerase (PAPD5 or PAPD7), and a substrate specific RNA-binding factor (ZCCHC7 or ZCCHC8). Several TRAMP-like complexes exist with specific compositions and are associated with nuclear, or nucleolar RNA exosomes. Nucleus, nucleolus. Event=Alternative splicing; Named isoforms=2; Name=1; IsoId=Q8N3Z6-1; Sequence=Displayed; Name=2; IsoId=Q8N3Z6-2; Sequence=VSP_013841, VSP_013842; Contains 4 CCHC-type zinc fingers. It is uncertain whether Met-1 or Met-2 is the initiator. Sequence=AAH36940.1; Type=Erroneous initiation; Sequence=BAB15418.1; Type=Misce

Physical Characteristics

Quantity: 1 ml (culture supernatant) or 100ug at 1mg/ml (purified)

Format: culture supernatant or purified material

Host/Isotype: mouse IgG1

Clonality: monoclonal; ID JH43.1.2F4

Formulation: culture supernatant contains 0.02% NaN₃. Purified material contains 30% glycerol, PBS and 0.02% NaN₃

Specificity: monospecific for human ZCCHC7; see "Microarray Analysis" below

Reactivity: human; not tested for cross reactivity in other species

Stability/Storage: 12 months long term: -20°C; short term: 4°C; avoid freeze-thaw cycles; aliquot as required

Handling Notes: small volumes of antibody may occasionally become entrapped in the seal of the product vial during shipment and storage; briefly centrifuge the vial on a tabletop centrifuge to dislodge any liquid in the container cap

Tested Research Applications

- **Western Immunoblotting - PASS by SOP**
- **Immunoprecipitation - Not recommended**

Antibody tested as purified IgG. Optimal dilution to be determined by user.

Quality Assurance

Notes:

1. Please refer to the SOP manual ([click here](#)) to for detailed explanation of the experimental setup and thresholds used to evaluate the results shown below.
2. Predicted MW of the target protein expressed in the western immunoblotting (WB) and immunoprecipitation (IP) experiment is **34.32kDa** and the three fusion tags (venus, 3xFLAG, V5) adds another **38.7kDa** to this protein. Therefore the protein is expected to migrate **~73.02kDa** in a denatured SDS-PAGE gel.
3. See results below for any applicable cautionary notes.

Guanajibo Research and Innovation Park
4005 St B Road 114 Km 1.3
Mayaguez, PR 00682

T 787.806-4100 Ext 233
F 787.806-4006
www.cdi-lab.com

© CDI Laboratories 2017.
For research use only. Not for use in humans or for diagnostics.
Lot-specific COA version tracker: v1.0.0



Anti-Human ZCCHC7, monoclonal

Alternate Names: Zinc finger CCHC domain-containing protein 7

Cat. No. 14-020

FOR RESEARCH USE ONLY


NOT FOR USE IN HUMANS



NEXTGEN PROTEOMICS

PRODUCT DESCRIPTION

Quality Assurance (HuProt™ Array)

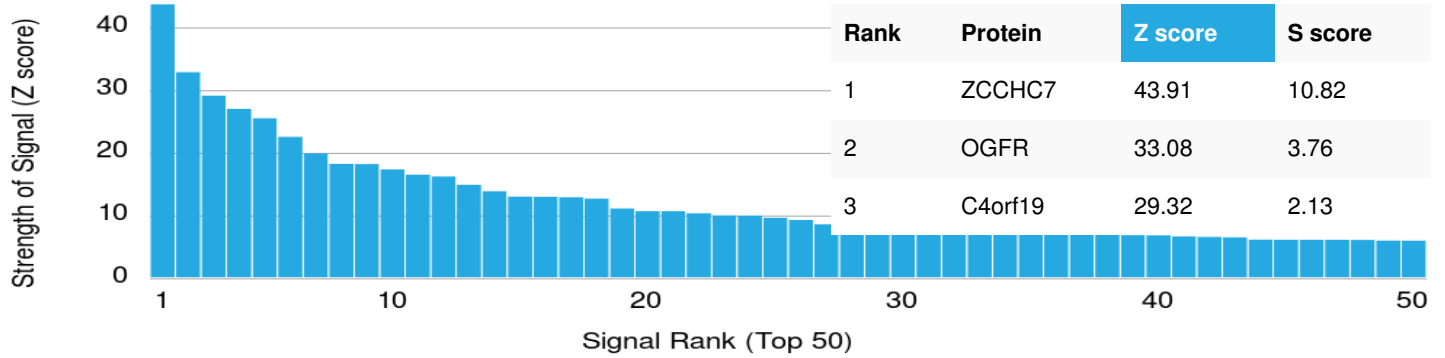

About Z and S Scores: The z-score represents the strength of a signal that a mAb (in combination with a fluorescently-tagged anti-IgG secondary antibody) produces when binding to a particular protein on the HuProt™ array. Z-scores are described in units of standard deviations (SDs) above the mean value of all signals generated on that array. To be considered a definite binding event between a mAb and a protein target, the z-score must have a value of at least 2.5 (SDs above the mean). The s-score represents the relative target specificity of a mAb, and is the difference (also in units of SDs) between the z-score a mAb generates upon binding to one protein and the next highest z-score that the mAb generates against another protein. For example, if a mAb binds to protein X with an z-score of 43 and to protein Y with an z-score of 14, then the s-score for the binding of that mAb to protein X is equal to 29.

Guanajibo Research and Innovation Park
4005 St B Road 114 Km 1.3
Mayaguez, PR 00682

T 787.806-4100 Ext 233
F 787.806-4006
www.cdi-lab.com

© CDI Laboratories 2017.
For research use only. Not for use in humans or for diagnostics.
Lot-specific COA version tracker: v1.0.0



NEXTGEN PROTEOMICS

Anti-Human ZCCHC7, monoclonal

Alternate Names: Zinc finger CCHC domain-containing protein 7

Cat. No. 14-020

FOR RESEARCH USE ONLY


NOT FOR USE IN HUMANS



PRODUCT DESCRIPTION

Quality Assurance (continued)

WESTERN BLOT - Predicted MW 73.02kDa

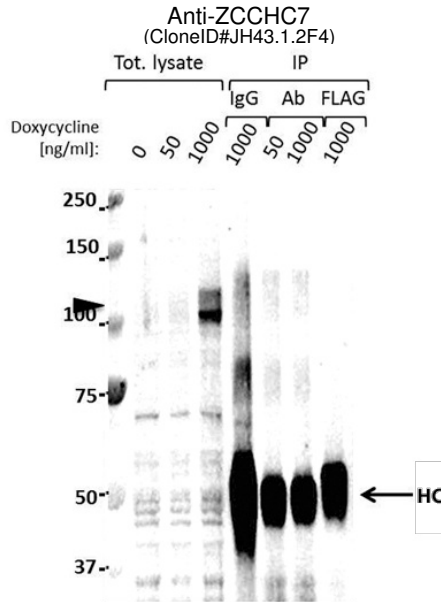

Tet-ON HeLa cells were transfected with construct encoding ZCCHC7 (BC034022.1) with an N-terminal fusion of FLAG, YFP (Venus) and V5 tags, under a tet-inducible promoter. These cells were stimulated with 0, 50 or 1000 ng/ml doxycycline. Immunoprecipitation (IP) was carried out using 5 μ g of either IgG, CDI mAb Anti-ZCCHC7 (cloneID# JH43.1.2F4) or 1 μ g of FLAG-M2 (Sigma). Immunoblotting was performed using 0.2 μ g/ml CDI mouse mAb Anti-ZCCHC7 (cloneID# JH43.1.2F4). HC=Heavy chain.

Cautionary notes for IP and WB experiments:

1. Target protein runs substantially heavier or lighter than expected

Guanajibo Research and Innovation Park
4005 St B Road 114 Km 1.3
Mayaguez, PR 00682

T 787.806-4100 Ext 233
F 787.806-4006
www.cdi-lab.com

© CDI Laboratories 2017.
For research use only. Not for use in humans or for diagnostics.
Lot-specific COA version tracker: v1.0.0

